#### **Product Datasheet**

# GluR2 (phospho-Ser880) antibody

Catalog No: #11292

Package Size: #11292-1 50ul #11292-2 100ul



Orders: order@signalwayantibody.com Support: tech@signalwayantibody.com

## Description


Dodonption	
Product Name	GluR2 (phospho-Ser880) antibody
Host Species	Rabbit
Clonality	Polyclonal
Purification	Antibodies were produced by immunizing rabbits with synthetic phosphopeptide and KLH conjugates.
	Antibodies were purified by affinity-chromatography using epitope-specific phosphopeptide. Non-phospho
	specific antibodies were removed by chromatogramphy using non-phosphopeptide.
Applications	WB
Species Reactivity	Hu Ms Rt
Specificity	The antibody detects endogenous level of Glutamate receptor 2 (Precursor) only when phosphorylated at
	serine 880.
mmunogen Type	Peptide-KLH
Immunogen Description	Peptide sequence around phosphorylation site of serine 880 (I-E-S(p)-V-K) derived from Human Glutamate
	receptor 2.
Conjugates	Unconjugated
Farget Name	GluR2
Modification	Phospho
Other Names	GRIA2; GluR-2; GluR-B; GluR-K2; AMPA 2
Accession No.	Swiss-Prot: P42262NCBI Protein: NP_000817.2
SDS-PAGE MW	100kd
Concentration	1.0mg/ml
Formulation	Supplied at 1.0mg/mL in phosphate buffered saline (without Mg2+ and Ca2+), pH 7.4, 150mM NaCl, 0.02%
	sodium azide and 50% glycerol.
Storage	Store at -20°C for long term preservation (recommended). Store at 4°C for short term use.

#### **Application Details**

Predicted MW: 100kd

Western blotting: 1:500~1:1000

#### **Images**

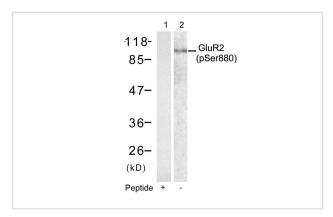

Western blot analysis of extracts from mouse brain tissue using Glutamate receptor 2(Precursor)(phospho-Ser880) antibody #11292(Lane 2) and the same antibody preincubated with blocking peptide(Lane1).

### Background

lonotropic glutamate receptor. L-glutamate acts as an excitatory neurotransmitter at many synapses in the central nervous system. Binding of the excitatory neurotransmitter L-glutamate induces a conformation change, leading to the opening of the cation channel, and thereby converts the chemical signal to an electrical impulse. The receptor then desensitizes rapidly and enters a transient inactive state, characterized by the presence of bound agonist.

Cull-Candy, S. et al. (2006) Curr. Opin. Neurobiol. 16, 288-297.

Hanley JG, et al. (2002) Neuron. 34(1): 53-67

Note: This product is for in vitro research use only and is not intended for use in humans or animals.